

## TITLE: Replacing the ME220 Quick Disconnect Insert

The Quick Disconnect Insert and Cap Assembly connects the Supply and Waste Bottles to the ME220 *Water Works* System.

### Included

Quick Disconnect Insert (2)

1. Disconnect the bottle by pressing the thumb latch on the Quick Disconnect Cap. See Figure 1.

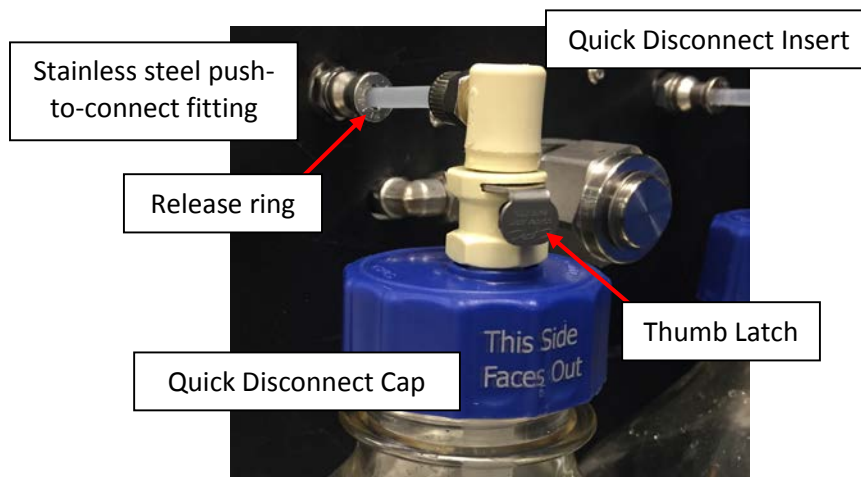


Figure 1 – Quick Disconnect Components

2. Disconnect the existing Quick Disconnect Insert by pushing the 'release ring' in, while pulling the tubing out. Discard the used insert.



Figure 2 – Removing the Quick Disconnect Insert

3. Connect the replacement Quick Disconnect Insert by pushing the tubing into the fitting, approximately 0.5", until you reach a hard stop.
4. Orient the bottle in the *Water Works* Bay with the label on the Quick Disconnect Cap facing out, as shown in Figure 3.



**Figure 3 – Connecting the Bottle to the *Water Works* System**

5. Reconnect the bottle by pressing the Quick Disconnect Insert into the coupling on the Quick Disconnect Cap until it clicks into place.